

# THE BATTERY-OPERATED CART WITH MECHANICAL HEIGHT ADJUSTER

pro-cart

Recorded information must be available everywhere at any time in order to ensure optimum and efficient care for patients. **pro-cart** has advanced to become the optimum assistant through the latest development in the form of an integrated power supply integrated within the base and independent of mains electricity.

Features of the supply / charging unit:

- Power supply of up to 3 items of equipment
- Fully enclosed and maintenance-free enclosure
- Battery with power of 12V and 55Ah
- Operating time of approx. 10 hours
- Input voltage: AC 115V/230V, 50Hz/60Hz
- Output voltage: DC 12V, 20V, 24V <sup>1)</sup>

Advantages of mechanical height adjustment:

- Variable adjustability (max. payload up to 16 kg)
- Available as a standing and seated workplace at any time and hence the guarantee for an ergonomically correct working position (proper attachment via mounted clamp levers)
- Control panel with additional information about charge status and operating voltage



The control panel is user friendly and attached at the front of the shelf to be easily visible. It provides information about the charging status.

<sup>1)</sup> Other output voltages are available upon request

# THE BATTERY-OPERATED CART WITH MECHANICAL HEIGHT ADJUSTER

**pro-cart**

Designation: **pro-cart** mobile computing with electrical height adjuster  
 Order number: **MC.2001.XXX** <sup>1)</sup>



Base 35 U (1120 mm), with 4 stabilizers, fitted with:

- 4 twin swivel castors with Ø 125 mm, with 2 lockable castors
- 1 mechanical height adjuster
- 1 shelf with control panel with display of charge status, dc / ac voltage and handles mounted at the side
- 1 monitor mount VESA 75 / 100 up to 14 kg monitor weight
- 1 keyboard drawer with pullout mouse pad right / left
- 1 push handle
- 1 battery 55Ah
- 1 supply / charging unit
- 1 connection cable for supply / charging unit <sup>2)</sup>

**Load capacity:** 80 kg / 176 lbs

**Finish:**

- RAL 7035 light grey power-coated
- Handle option see <sup>1)</sup>

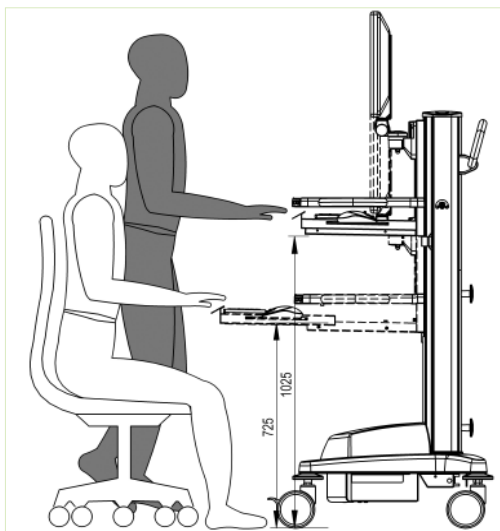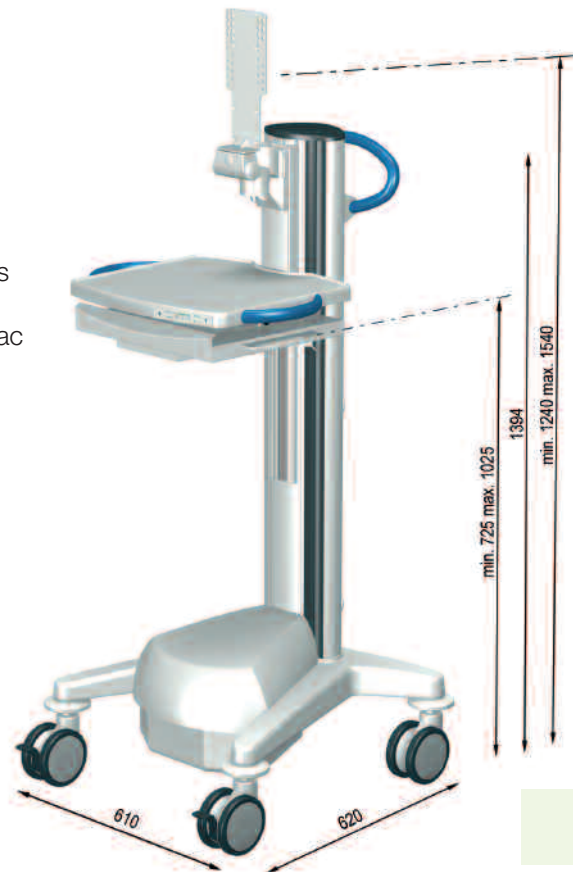

**pro-cart** can be adjusted to suit the appropriate working position. This feature ensures correct **ergonomic working conditions** at all times.



The monitor mount permits a tilt of  $\pm 15^\circ$  and can also be swivelled through  $\pm 20^\circ$ , and is suitable for monitors up to 14 kg with VESA 75 / 100 standard.

<sup>1)</sup> **Trim colour:** Make sure you replace .XXX by .903 = RAL 7035 light grey, .908 = RAL 5015 sky blue or .909 = RAL 1003 signal yellow  
<sup>2)</sup> The connection cable for connecting the supply / charging unit with the equipment must be ordered separately (ZV.9374.999)