




Medical cart
classic-cart

NEW PRODUCT

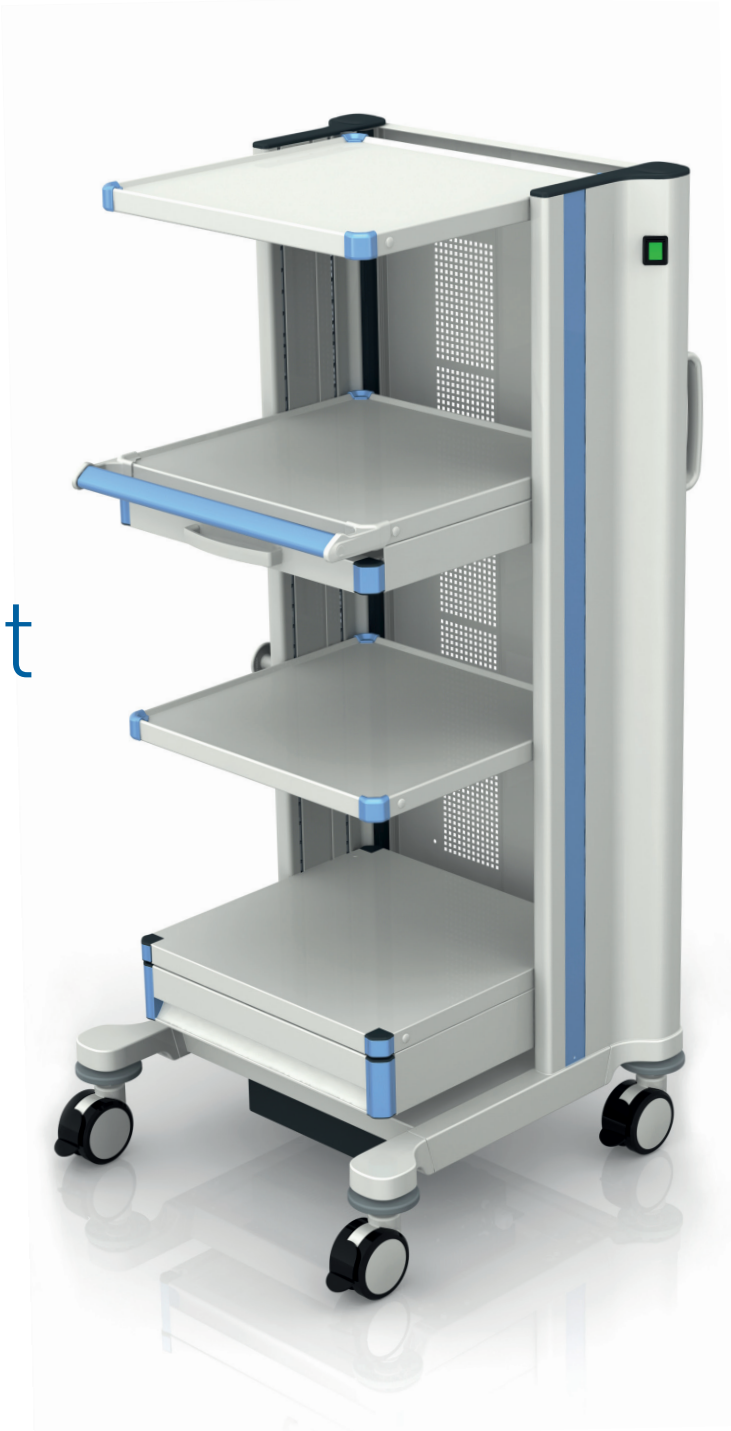

TECHNICAL DETAILS

Load capacity

Total load capacity	150 kg / 330 lbs
Load capacity of shelf	50 kg / 110 lbs distributed load

Dimensions - Basic frame

Width in mm	663 / 763 / 1263
Depth in mm	660
Total height in mm	840 / 968 / 1128 / 1288 / 1448 / 1608 / 1768

Dimensions - Shelf

Width in mm	500 / 600 / 1100
Depth in mm	491

Electrification

Isolating transformer	1200 VA / 1600 VA / 2000 VA
Central on / off switch	yes
Socket strips	yes

Get in touch »

Europe: +49 89 614425-0 | sales@itd-cart.com
America: +1 800 947-3901 | salesusa@itd-cart.com
China: +86 21 5031 8818 | saleschina@itd-cart.com
India: +91 86 0015 4064 | salesindia@itd-cart.com
Australia: +61 3 9330 5050 | salesaustralia@itd-cart.com


www.itd-cart.com



CUSTOMIZED DESIGN

New colours

BASE COLOR




RAL 9016
Traffic white

DECORATIVE TRIM COLORS



RAL 7016
Anthracite grey



RAL 9016
Traffic white

New design

PERSONALIZED CORPORATE DESIGN





DS.0270.799 -02/2026. Printed in Germany. All rights to make changes reserved.

HIGH QUALITY DESIGN WITH OPTIMIZED FUNCTIONALITY



New corners –
elegant design with improved
hygiene properties

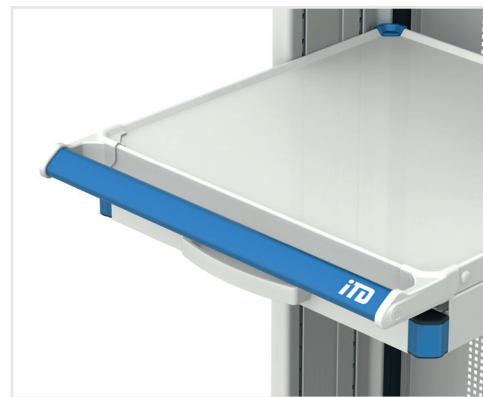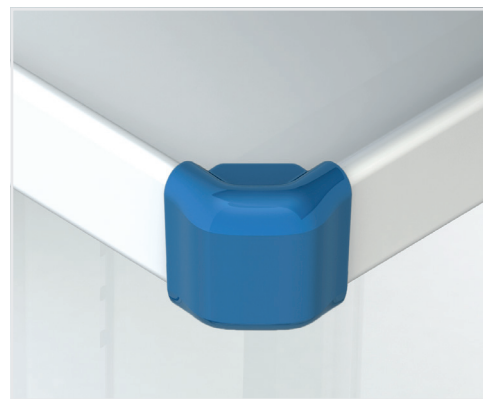

Push handle –
comfortable feel, customizable
to corporate brand standards

Single, smooth column –
offers excellent hygiene
features and a distinct design





Customized design –
highly recognizable

New arm position –
large range of motion



**Additional
rear handles –**
easy align-
ment of the
cart from any
position

International socket strips –
unrestricted accessibility

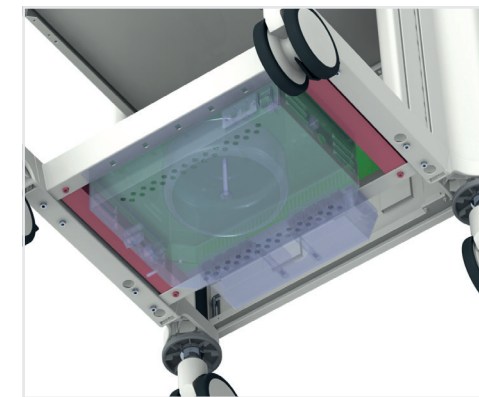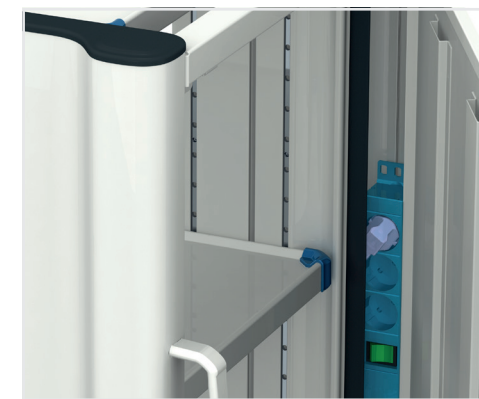

Counterweights –
invisibly integrated and flexibly
combined with isolating transformers

Air-permeable rear door –
prevents heat congestion of the
equipment and unwanted access

